

SEESHEEN



SOLE DISTRIBUTOR



BUILT FOR PORTABILITY



New Generation
Image Processor + Light
Source+Monitor


Image Processor



VIDEO ENDOSCOPE IMAGE PROCESSOR HV-3101


10 INCHES

Portable Endoscope System

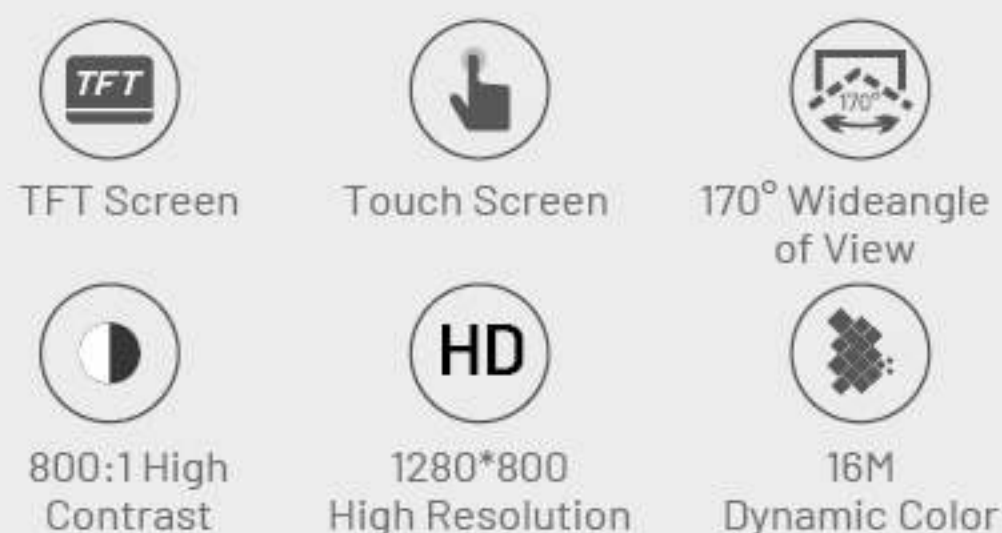


EXCELLENT CHOICE OF BEDSIDES USE

Better user interface

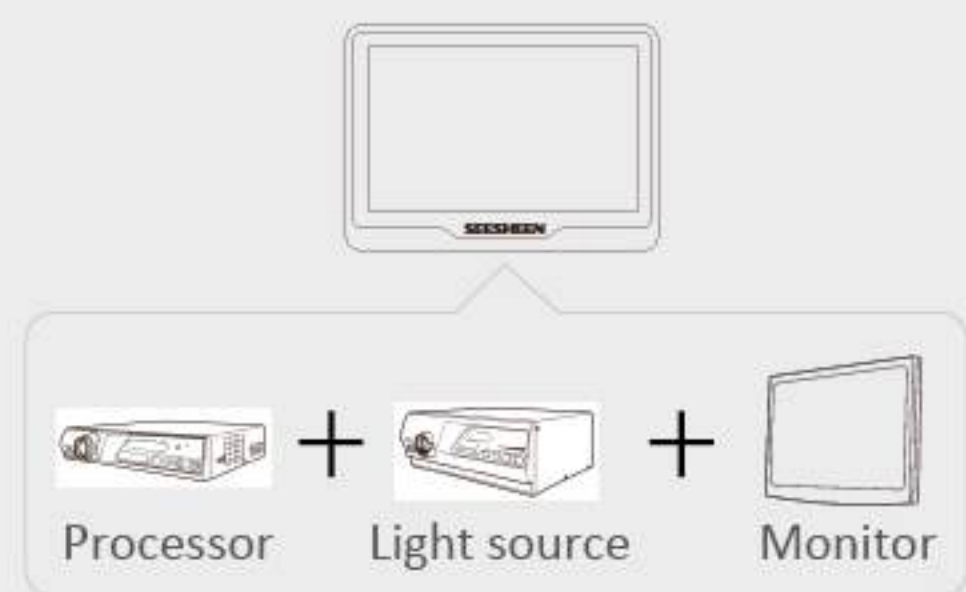

No delay feeling of touch screen operating
Allow operating with gloves

All in one Highly Integrated Design

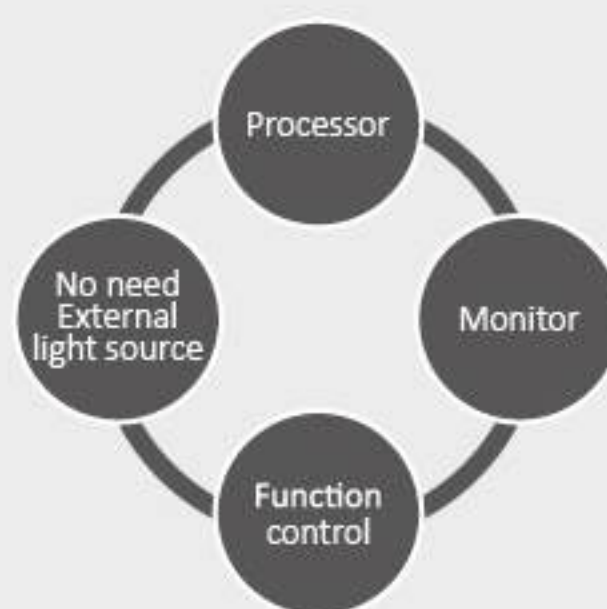

Detachable rechargeable battery supporting 4 hours working



Intelligent split screen display more efficient



Detachable SD card



HDMI+BNC (CVBS) output

CMOS+LED

FLEXIBLE VIDEO
CHOLEDOCHOSCOPE

SUCTION BUTTON FIXTURE

Firmly fasten the suction button



PROTABLE

Compact design for less
space in O.R.

CMOS + LED TECHNOLOGY

Non-fiber design makes it durable and
easy for maintenance

Mainly applied for diagnosis and treatment of gallstones, such as PTCS surgery. The flexible insertion tube can pass a small nvasive puncture hole to the common bile duct and hepatobiliary.



Model	Insertion Tube	Channel	Working Length	Bending (Up/Down)	Field of View	Depth of Field
CH-1349	4.8 mm	2.0 mm	380 mm	160°/130°	120°	3-50 mm

SEESHEEN




SOLE DISTRIBUTOR



UR-1328 SERIES
SINGLE-USE FLEXIBLE
CHOLEDOCHOSCOPE

Dual light sources, dedicated irrigation, large 2.0mm working channel.

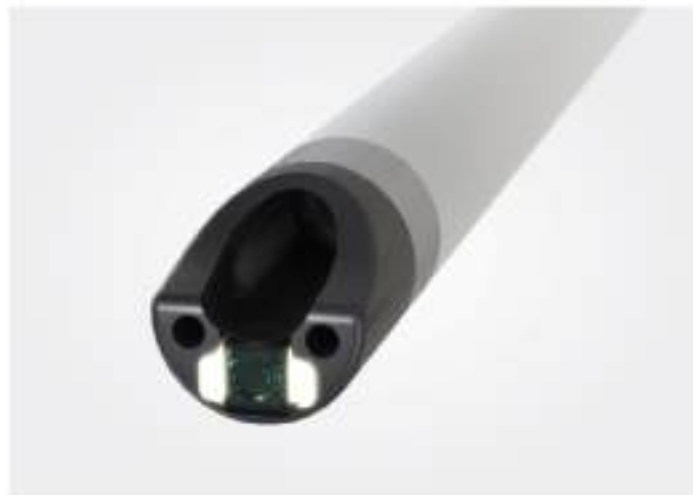

Ergonomic and lightweight, with deflection lock reduces fatigue.



Provides four-way steering, provides traditional sens of operation.



Attaches to duodenoscope for single operator use.



EXCLUSIVE CHANNELS FOR DIFFERENT OPERATIONS



Instrument Irrigation flow control



Suction controlled



Suction



Water jet/ Guide wire insertion



Reduces cost and risk associated with reprocessing and cross-contamination

Product Parameter

Model	Insertion Tube	Working Channel	Water jet Channel	Suction Channel	Working Length	Bending (U/D/L/R)	Field of View	Depth of Field
CHD-1	3.7 mm	2.0 mm	0.6 mm	0.6 mm	2120 mm	120°	120°	3-50 mm

Image Processor

HV-3101
Image Processor

10.1 inches Touch screen

Detachable rechargeable battery supporting 4 hours working


Detachable SD card

Portable and lightweight design

HD + AV output

HD IMAGE

FLEXIBLE VIDEO CYSTONEPHROSCOPE

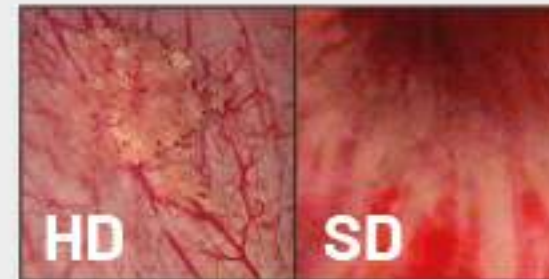

HD IMAGE

HD image gives more details of bladder

CONTROL BUTTON

1. Zoom
2. Capture/Recording
3. Freeze



CMOS + LED TECHNOLOGY

Fiber-free design makes it more durable



Insertion tube: 15.6Fr

Working channel: 6.6Fr

PORTABLE

Equipped with image processor HV-3101
Support external monitor



APPLICATION

Bladder examination, therapy and kidney calyx observation during PCNL.

Model	Insertion Tube	Channel	Working Length	Bending (Up/Down)	Field of View	Depth of Field
CY-1356	15.6Fr	6.6Fr	380 mm	210°/130°	90°	3-50 mm

UR-1328

FLEXIBLE VIDEO URETERORENOSCOPY



EXCELLENT MANEUVERABILITY

1. Zoom
2. Capture/Recording
3. Freeze

4-way angulation for accurate positioning



CMOS + LED TECHNOLOGY

Fiber-free design makes it more durable



Insertion tube: 7.5Fr

Working channel: 3.6Fr

Access sheath: ≥ 9.5 Fr

PORTABLE

Equipped with image processor HV-3101
Support external monitor



APPLICATION

Bright and clear image for non-invasive kidney lithotripsy: retrograde intrarenal surgery.

RIRS

Model	Insertion Tube	Channel	Working Length	Bending (Up/Down)	Rotation (Left/Right)	Field of View	Depth of Field
UR-1328	7.5Fr	3.6Fr	670 mm	275°/275°	120°/120°	120°	3-50 mm

DISPOSABLE

FLEXIBLE VIDEO URETEROSCOPE



INTRODUCTION

Single use flexible video ureterorenoscope delivers high quality image for visual inspection with a clean and safe clinical flow. It's affordable without reprocessing and maintenance, and it benefits both patient and hospital by effectively avoiding cross contamination.

PRODUCT FEATURES

- | Excellent image quality with highly integrated CMOS and LED
- | 4-way bending facilitates multidimensional observation
- | Angle-lock function by resistance force gives precise angle control and reduces operating fatigue
- | 3.6 Fr working channel ensures good perfusion flow and reduces the risk of surgical infection



Capture/Recording
Freeze



Compound Bending
Up/Down 275°



Slim Diameter



Left/Right Rotation



Angle-Lock Switch